

**Quadruple layer program to boost skin care potential**

# SKIN CURLING



WITHBEAUTY



**KC Certified**

KC 인증서 보유



**CE Certified**

CE 인증서 보유



**Patent Application**

제10-1772689호

특허 열을 이용한 미용관리기



**Patent Application**

제10-1763927호

열을 이용한 미용기기용 헤드

## **Skin Curling's Skin quadruple strengthen program Quadruple Program**

How to braking a skin watch that can not be disobeyed ?  
Skin Curling's quadruple program changes the depth of skin care and helps to make healthy and resilient skin from face to body.





Program 1

## Heating

The warmth of 42 °C that erases the memory of aging through wrinkle improvement of skin



Program 2

## Electro therapy

Microcurrents to enhance whitening, anti-aging and lifting effects



Program 3

## Vibration

Power motion that pat skin 6000 times per minute



Program 4

## 3-Color therapy

Three kinds of color therapy that can be selected to fit the skin condition



Program 01

## Heating

The warmth of 42 °C that erases the memory of aging through wrinkle improvement of skin


The 42 degree heat promotes the production of **HSP**, helping to reduce strikingly wrinkles and regenerate skin cells which are damaged by stress and ultraviolet rays. In addition, a warm feeling like a steam towel promotes **blood circulation in the skin**, helping to make the ideal skin for skin care.

### \* HSP(Heat Shock Protein)

It is a Heat Shock Proteins that contributes to restoring damaged skin due to stress and external environment and activating immune cells.

In particular, when the skin temperature is raised to 42 degrees, the HSP in the cell is increased.

HSP is effective in inhibiting wrinkles by reducing the expression of collagenase.

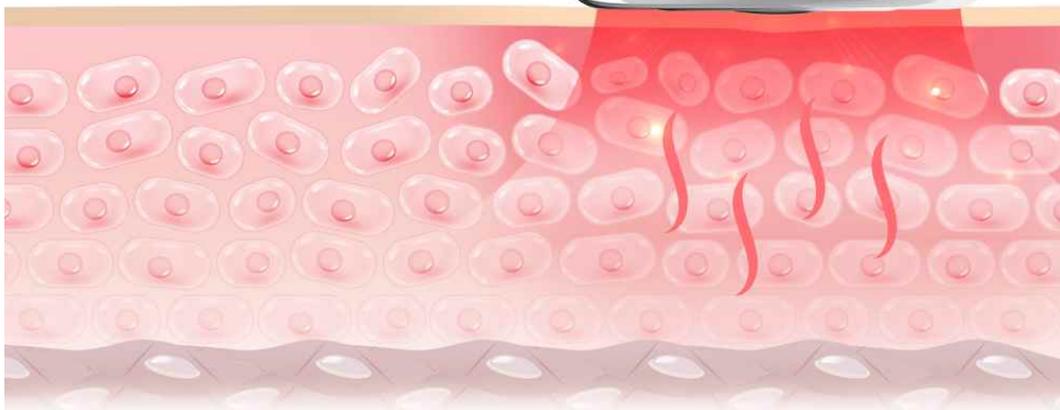

Program 02


# Electro Care

**Microcurrents** to enhance whitening, anti-aging and lifting effects


Electro Care is a method of **delivering cosmetic active ingredients deep into the skin** without irritation through electricity.

When electricity is transferred to the skin, spacing occurs between the skin cell tissues, allowing active ingredients penetration and waste discharge according to the flow of the electrons.


Skin Curling helps to reprogram the future of the skin by adding **implementing an excellent waveform** for wrinkle improvement and elasticity enhancement in the conventional ion introduction method that allows the cosmetic active ingredient to reach the dermis layer.

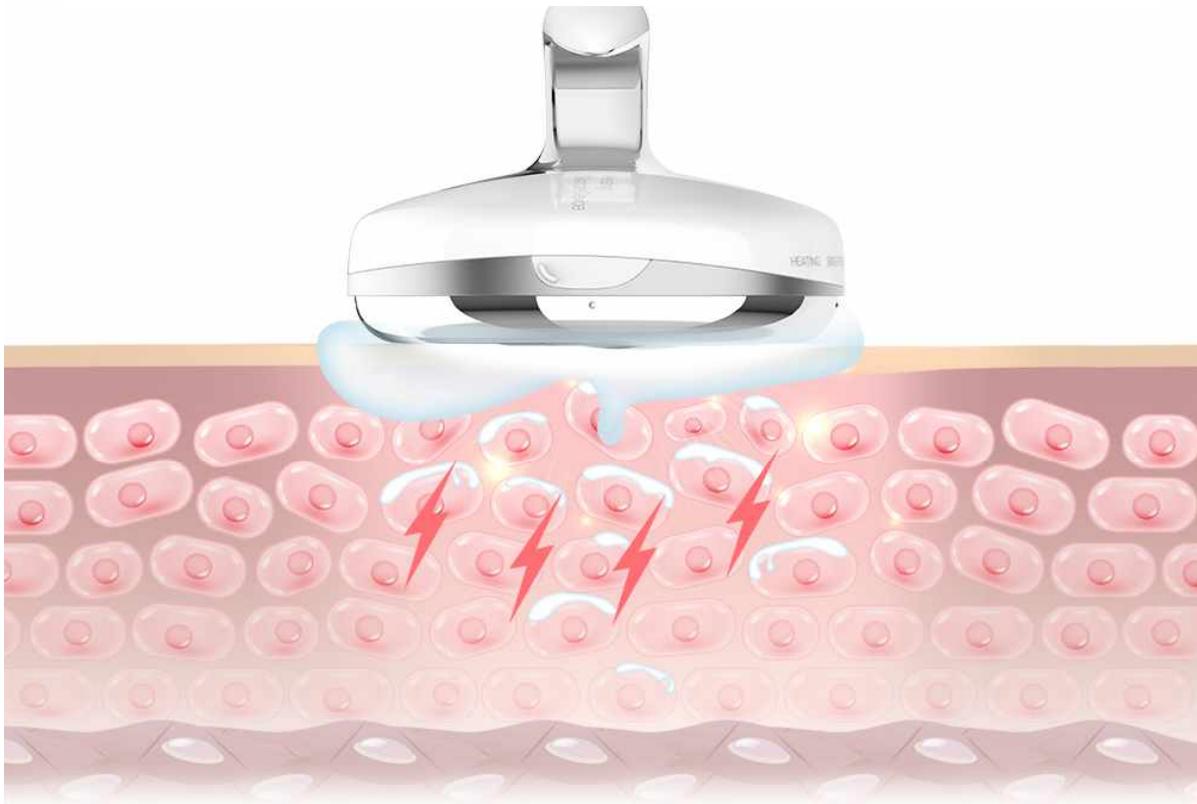


< Microcurrent waveform by mode >

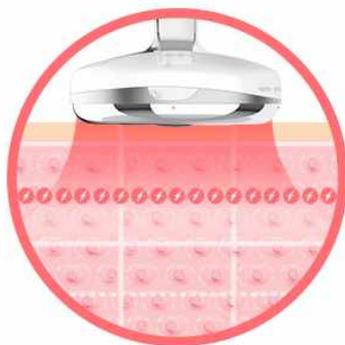

Brightening mode

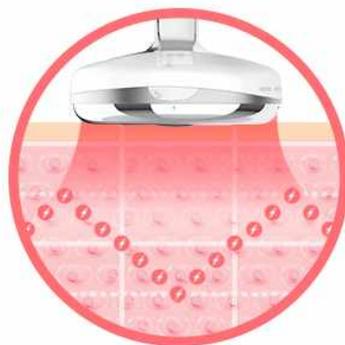

Anti-aging mode

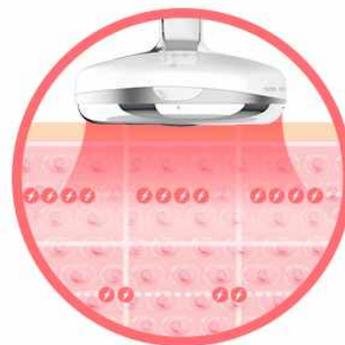


Lifting mode

Program 03

# Vibration

**Power motion** that pat skin 6000 times per minute

Skin Curling's power motion stimulates the capillaries and muscle layers of the skin and helps **improve blood circulation**. Power massage, which vibrates more than 6,000 times per minute, pat the skin so that the active ingredients of cosmetics that are not easily absorbed by hand can reach the dermis layer.

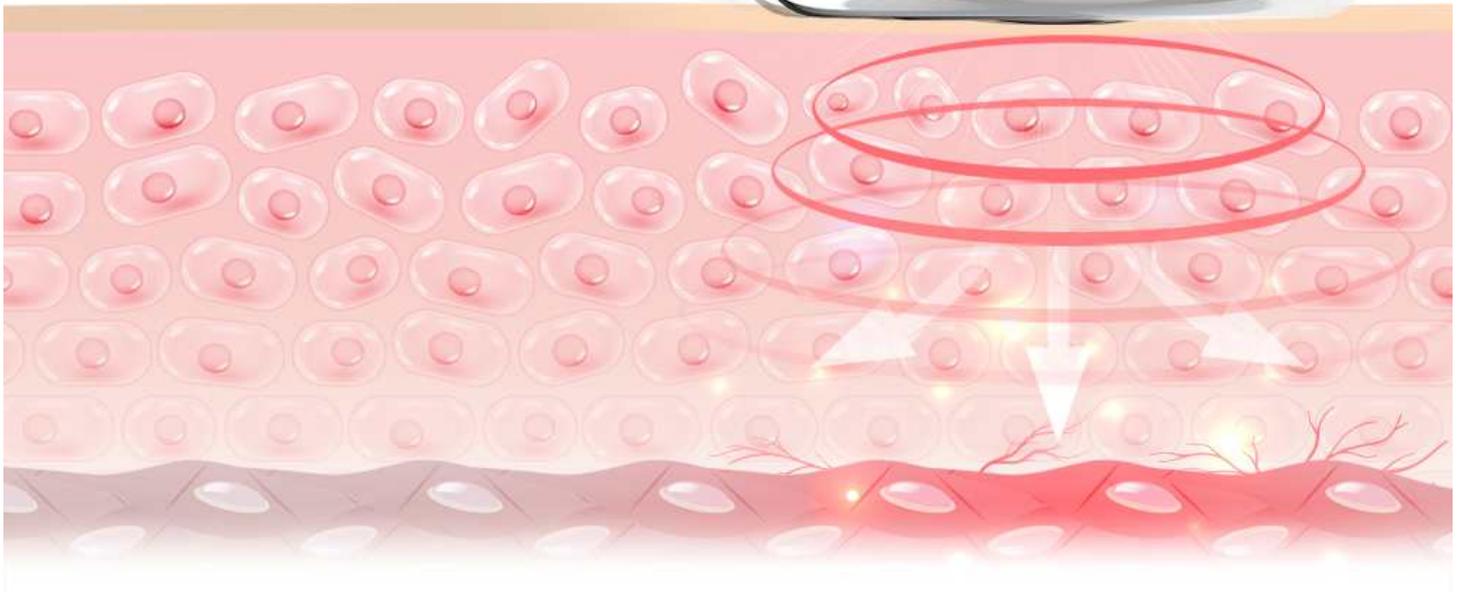

Program 04

# Color LED

**Three kinds of color therapy** that can be selected to fit the skin condition

**Aesthetic LED Masks have blind spot**

Skin Curling can gently roll the entire face and shine the light directly, so it can manage closely the major trouble spots for the skin without missing spots.

**Three wavelengths of light** that improve the skin can help whitening, wrinkle improvement, and elasticity to create clear and lively skin.



Color Led Mode 1

## Brightening Mode (GREEN)

Green Therapy, which helps reduce the pigmentation and blemish of the skin and helps to brighten the skin

Color Led Mode 2

## Anti-Aging Mode (RED)

Red Therapy, which promotes collagen production and helps improve skin elasticity




Color Led Mode 3

## Lifting Mode (BLUE & RED)

Red & Blue Double Action Therapy to help strengthen skin elasticity and calm skin

## <The Effect of three lights on each wavelength>

Depending on the wavelength of light, the degree of absorption to the skin is different and the effect also is different.

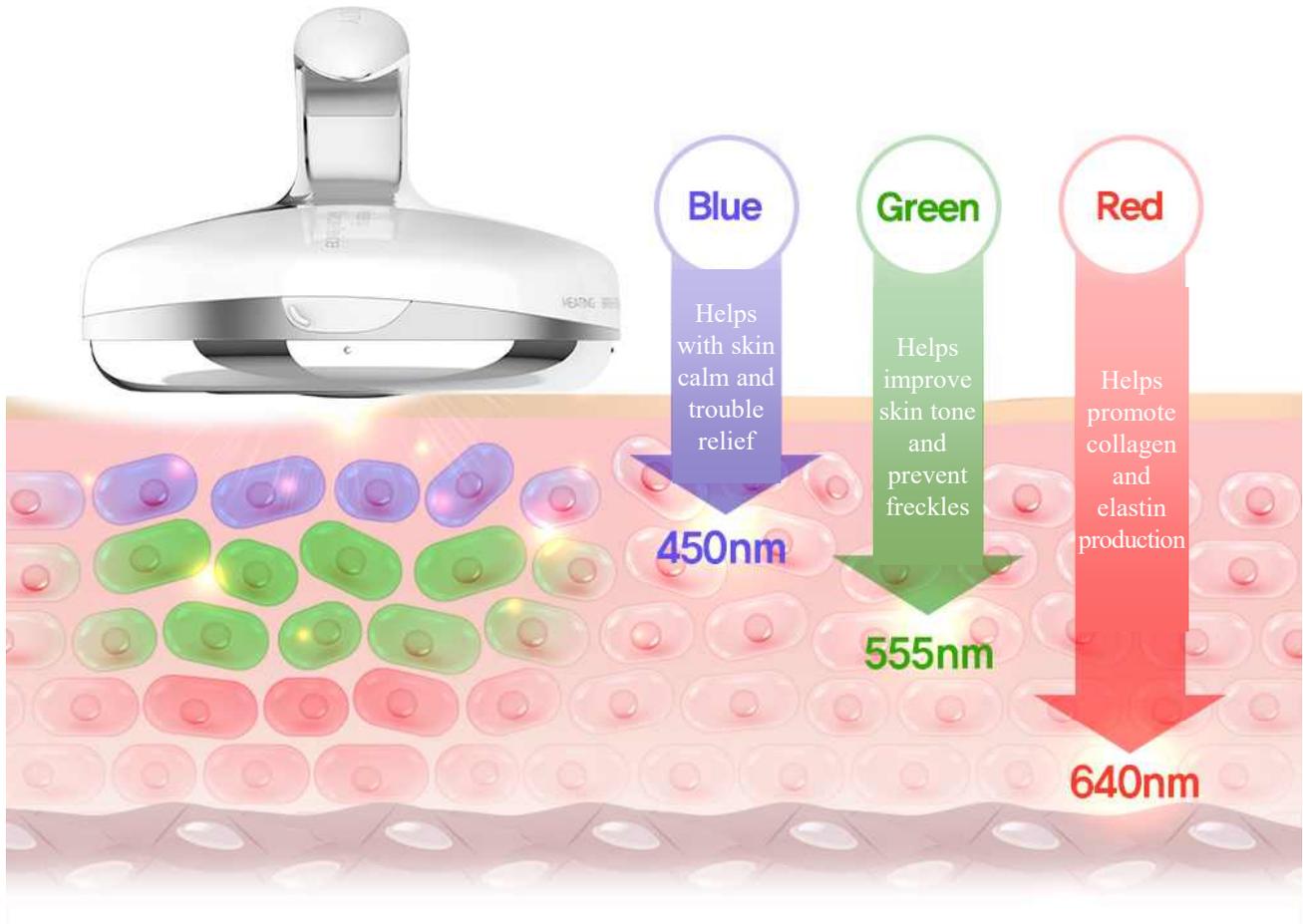

## <The Effect of Light on Skin>



Light energy infiltrated into the skin absorbs into the cells in the skin

Cell activation by light energy

Helps the biological energy metabolism in cells to be active

Light penetrated into the skin promotes cell activity and **helps normalize cell function** such as **skin regeneration** and **blood circulation improvement**, thereby helping to activate biological energy metabolism

## The cool waving head that is in full contact with every bend in the face and body

The triple waving head structure in full contact with the face and body, and gently massages muscles and blood vessels.

It is economical because cosmetics do not get pushed out of the skin even when it is rolling for a long time through Lock-up system that coop up in the spread massage cream on the face into the curved surface of the head.




## Rhodium plating that can be used safely for any skin type

Consumers' selection criteria for products that reach directly to the skin have become more complicated like that they consider whether it causes allergies to the skin or not harmful to the skin.

Skin Curling uses a brilliant 'Rhodium Plating' on the head area where the skin touches, satisfying the unique luxury and stability at the same time.

**Rhodium:** It is the platinoid elements harmless to the human body and is used for expensive luxury jewelry and watches  
**Tritan(BPA FREE):** It is an eco-friendly material without bisphenol, an environmental hormone, and is used in baby bottles and water bottles

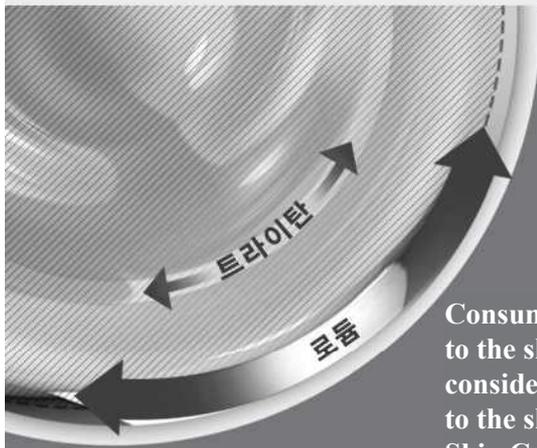

## Triple temperature control system for safe use

Skin Curling has applied a triple temperature control device to its product because we think safety is our top priority.

### **Thermistor**


Maintain 42 degrees by resisting and temperature error control


### **Bimetal**

Automatically shut down power when the temperature inside the product rises

### **Automatic control circuit**

After switching to standby mode, the operation of the product is automatically stopped when there is no mode change.

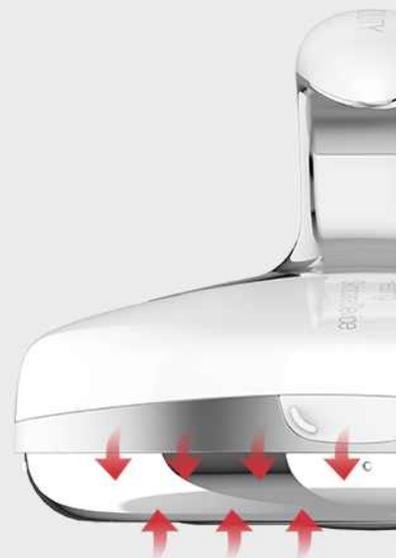

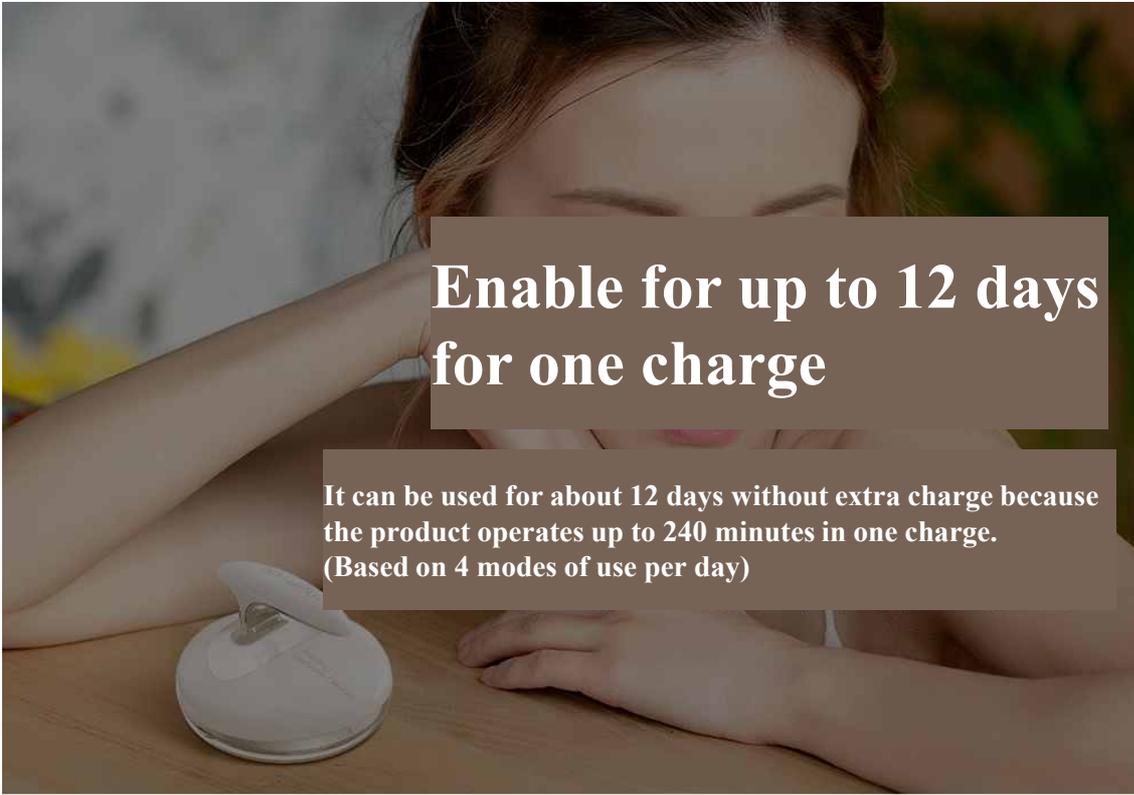

**Enable for up to 12 days  
for one charge**

It can be used for about 12 days without extra charge because the product operates up to 240 minutes in one charge.  
(Based on 4 modes of use per day)

## **Total Care System**

Total care solutions that can be managed from face to body at once

**It is designed to pass smoothly along the face and body lines,  
like a stone in a Curling game that passes over the ice.**



# Body Care TIP

This helps to create a resilient and neat body by smoothly making the bumpy line by rolling massage with body slimming cream.

**Back, Forearm, Abdomen, Calf, Butt**




**Dedicated  
pouches for  
portability and  
safety**

Dedicated pouches enhance product portability and ensure safe storage of products from external impacts and scratches.

# How To Use

It is easy to use because it is a button operation method.

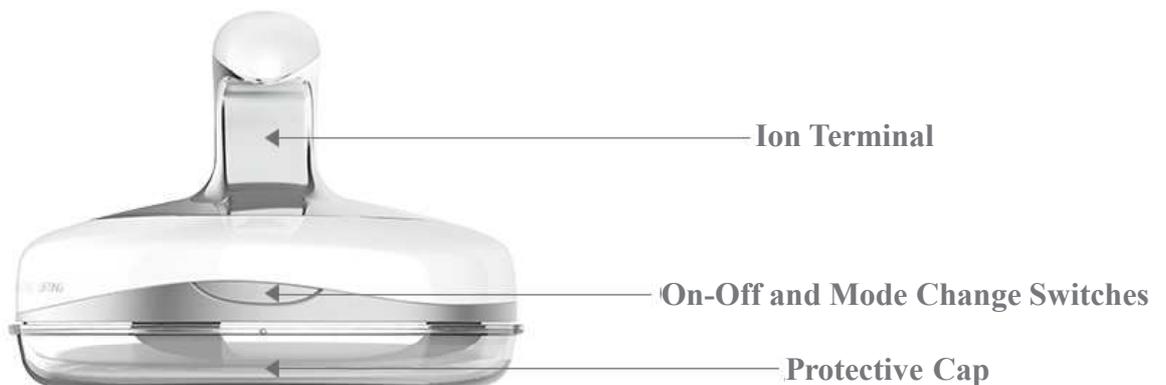

- 01** For galvanic effect, it use by your finger against to the ion terminal located at the top.
- 02** Press and hold the button until you hear a beep, then turn on the power and select the desired mode. Each time the button is pressed, the color of the lamp changes in the order of [ White ► Green ► Red ► Red & Blue ]
- 03** You put Skin Curling against to your skin and then massage it as if it were sliding.
- 04** If you want to stop using it during a massage, press and hold the button to turn off the power.
- 05** Clean the remaining cosmetic residue on the head, cover the head with a protective cap, and store.

**- TIP: If you use this product in state of the mask pack on your face, it can increase the absorption rate of the active ingredient in cosmetics.**

# Usage by mode

Renew your skin every day! Choose a mode that suits your skin troubles to complete your bright, elastic skin



## ☉ Heating

- 01 When you press the switch 1, the white light lights up with the "beep" sound.
- 2 Massage your arranged skin evenly with heated equipment.

**Total time: 3 minutes and 30 seconds**



## ☀ Brightening

- 01 When the switch is pressed twice, the green light is turned on.
- 2 Massage all over the skin like a circle after applying whitening cream to the skin

**Total time: 5 minutes**



## ☀ Anti-aging

- Thermotherapy + Microcurrent + Vibrating + Color LED (Red)
- 01 When the switch is pressed 3 times, the red light is turned on.
  - 2 Massage the area evenly after applying wrinkle-improving cream on forehead or eye area.

**Total time: 4 minutes and 30 seconds**



## ☀☀ Lifting

- 01 When the switch is pressed 4 times, the red and the blue lights are alternately turned on.
- 2 Massage the area evenly after apply elastic cream to the jawline or clown area, etc.

**Total time: 5 minutes**

# Names and components of each part

## //////////////////////////////////// The Names of each part //////////////////////////////////////



## //////////////////////////////////// Components //////////////////////////////////////



# Product information

## WITHBEAUTY PROUDCT SPEC

<b>Name</b>	Skin Curling / WB-HD-413A
<b>Product Composition</b>	Main Body, 1 Rechargeable Adapter, Manual
<b>Power Consumption</b>	8W
<b>Size</b>	72 x 72 x 52 (mm)
<b>Weight</b>	100g
<b>Materials</b>	A,B,S, Stainless Steel, Cu, Cr3+, Titanium & Rhodium Coating
<b>Charge Power</b>	DC 5V / 1,2A
<b>Power Supply</b>	Li-polymer 3.7V / 900mAh
<b>Certification No</b>	R-REM-WB9-WB-HD-411A
<b>Manufacturer</b>	Withbeauty Co., Ltd.
<b>Manufacturing Country</b>	Korea
<b>A/S</b>	(주)위드뷰티 1661-1981